


IS-III active S-narrow

2019.04 E Ver.01





Find Neobiotech online

Super Narrow
for Anterior Mandible
(Narrow Ridge)

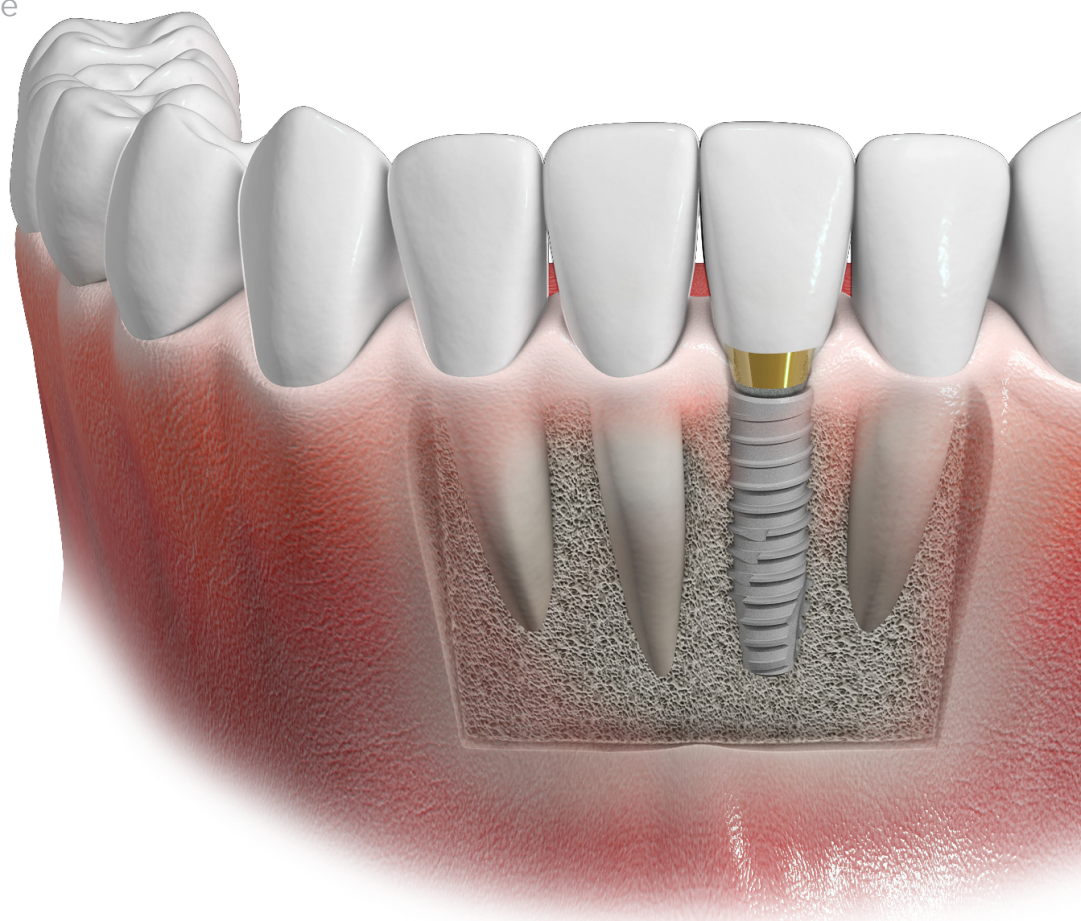


IS-III active S-Narrow Features

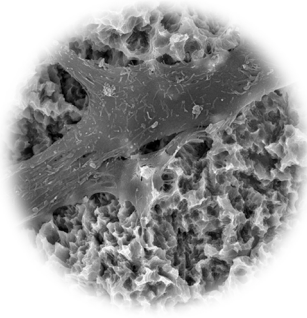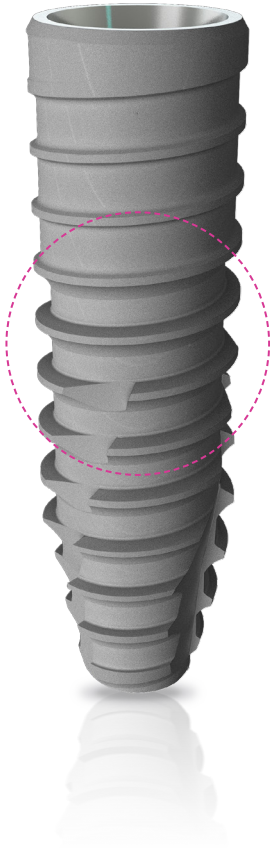


S.L.A. Surface

The new S.L.A Surface with 40% greater surface area and 50% more cell adhesion promotes osseointegration

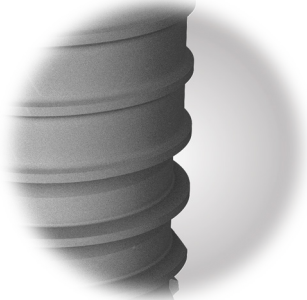


Straight Body

Easy to adjust implantation depth

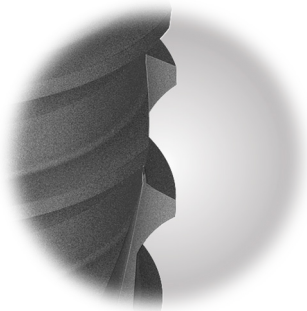

Deep & Wide Pitch

Optimum Pitch for osseointegration (0.9Pitch)

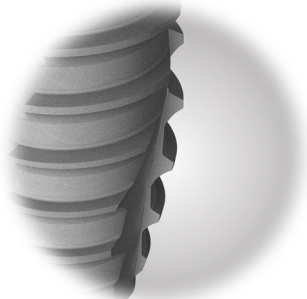

Wide Cutting Edge

Wide cutting edge and enlarged surface area enhances initial fixation and offers clinicians more stable implant placement


M2.0 Connection



One Piece Type

- Cover Screw
- Healing Abutment

Internal Conical Seal 11° 2.1 Hex

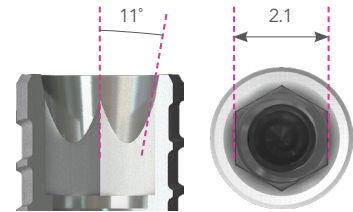


M1.6 Connection

Two Piece Type

- PickCap Impression Coping
- Temporary Abutment
- Cement Abutment
- Angled Abutment




Grade23 : Ti-6Al-4V ELI

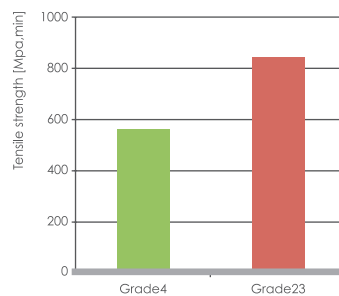

IS-III active S-Narrow Line Up

Fixture

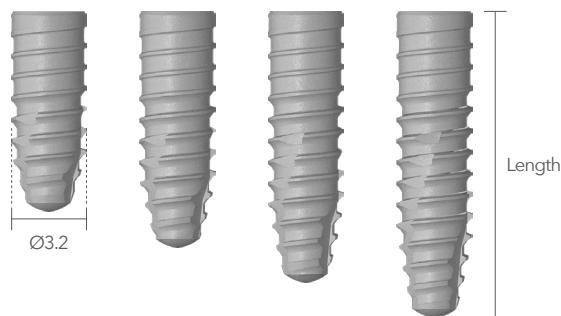

Diameter	Length	Code
Ø3.2	8.5 mm	IS3SN3008AP
	10.0 mm	IS3SN3010AP
	11.5 mm	IS3SN3011AP
	13.0 mm	IS3SN3013AP

* Cover Screw not included


Cover Screw & Healing Abutment

Cover Screw



Diameter	Code
Ø3.0	ISC25S

Healing Abutment

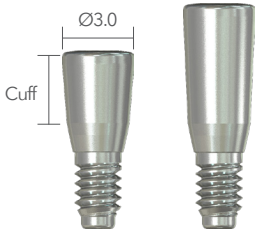

Diameter	Ø3.0	
Cuff	3 mm	5 mm
Code	ISH3003	ISH3005



Diameter	Ø3.8	
Cuff	3 mm	5 mm
Code	ISH3803	ISH3805

Prosthetic Flow Chart



Temporary
Abutment



Cement
Abutment

Custom
Abutment



Lab Analog



PickCap Impression
Copping



IS-III active
S-narrow (ø3.2)

IS-III active S-Narrow Prosthesis

Cement Abutment (Hex)




Diameter	Post	Cuff	Code
Ø3.5	7 mm	3 mm	ISAH3537
		4 mm	ISAH3547
		5 mm	ISAH3557
		7 mm	ISAH3577

Cement Abutment (Non- Hex)




Diameter	Post	Cuff	Code
Ø3.5	7 mm	3 mm	ISAN3537
		4 mm	ISAN3547
		5 mm	ISAN3557
		7 mm	ISAN3577

Temporary Abutment

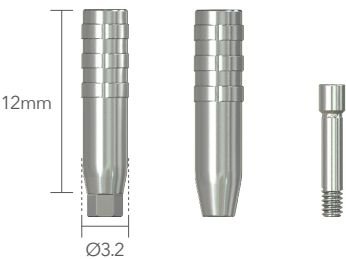

Diameter	Length	Connection	Code
Ø3.2	12 mm	Hex	ISAHT300
		Non Hex	ISANT300

PickCap Impression Coping

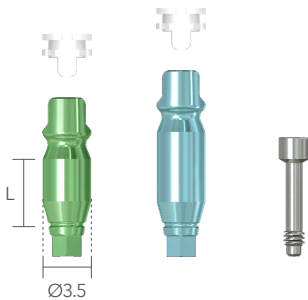

Diameter	Length	Code
Ø3.5	5 mm	ISSNPCIB3505C
	7 mm	ISSNPCIB3507C

* Screw & Cap included

Lab Anlog



Diameter	Code
Ø3.2	ISLA300


Surgical Kit

The top, middle, bottom tray and components of the S-narrow Kit is sold separately


Kit Tray Top	CSTKITTRAY1
Kit Tray Middle	SNKITTRAY2
Kit Tray Bottom	CSTKITTRAY3




Drill Extension

	Code	
	DE01	


Cortical Tap

	Diameter	Code
	Ø3.2	CT32S09


Hex Driver

	Code	
	HD1215S	


Guide Point Drill

	Diameter	Code
	Ø1.35	PD13


Cortical Drill

	Diameter	Code
	Ø3.1	ISCD32


Direction Pin

	Diameter	Code
	Ø3.6	DP36


Surgical Drill

	Diameter	Code
	Ø2.2	TSDS22CL
	Ø2.7	TSD27CL


Fixture Driver

	Code	
	ISFDSNRL	
	ISFDSNCL	


Connector

	Code	
	RC15S	

Drill Stopper


			
8.5mm	10mm	11.5mm	13mm
DS085CS	DS100CS	DS115CS	DS130CS

Torque Ratchet

	Code	
	TW60	

Surgical Kit

Tools from S-narrow Kit can be mounted on spare parts of the Neo Master Kit (Sold separately)

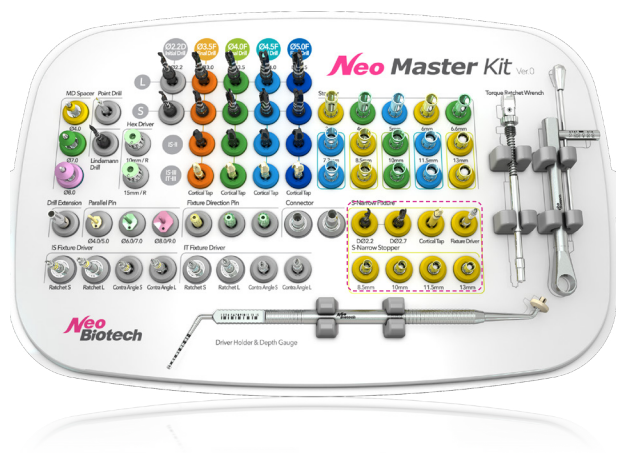

Surgical Drill



Drill Stopper



Cortical Tap




Fixture Driver



Drill Length

Surgical Drill

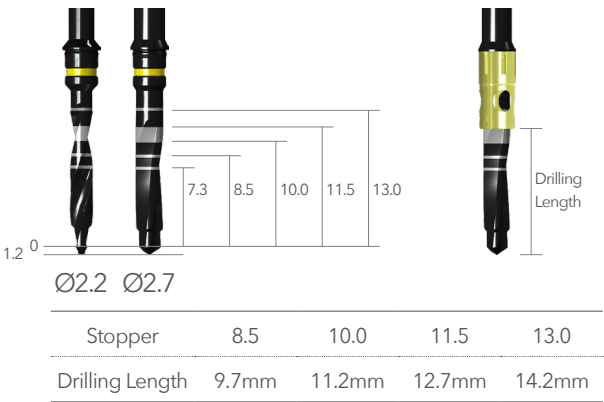

Cortical Drill




Drilling Protocol


Drilling Speed & Torque

Guide Point Drill, Surgical Drill : 1,200rpm / 35~45Ncm Cortical Tap : 50rpm / 50Ncm

IS-III active Fixture $\text{Ø}3.2 \times 10\text{mm}$ (D1/D2 bone)

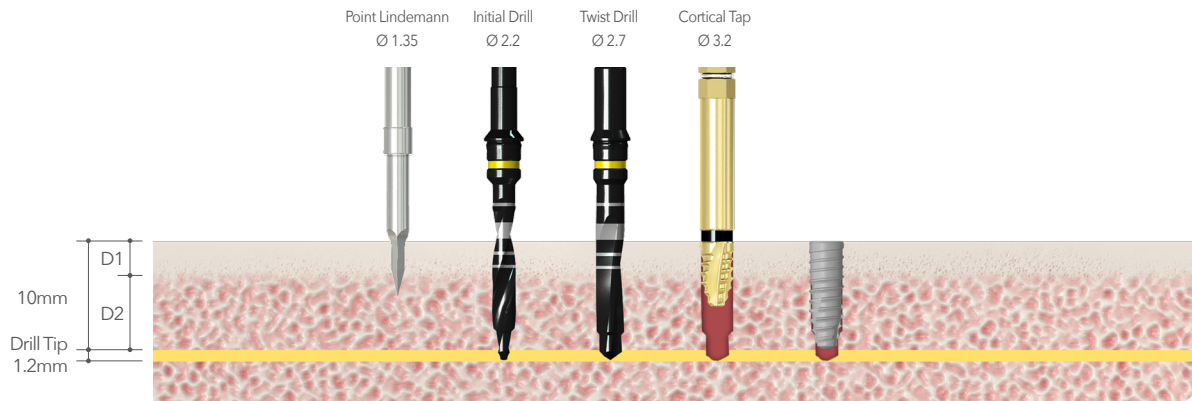

In case of soft(D4) bone, or in condition of getting fixed at implant apex area, $\text{Ø}2.2$ surgical drill is the final drill

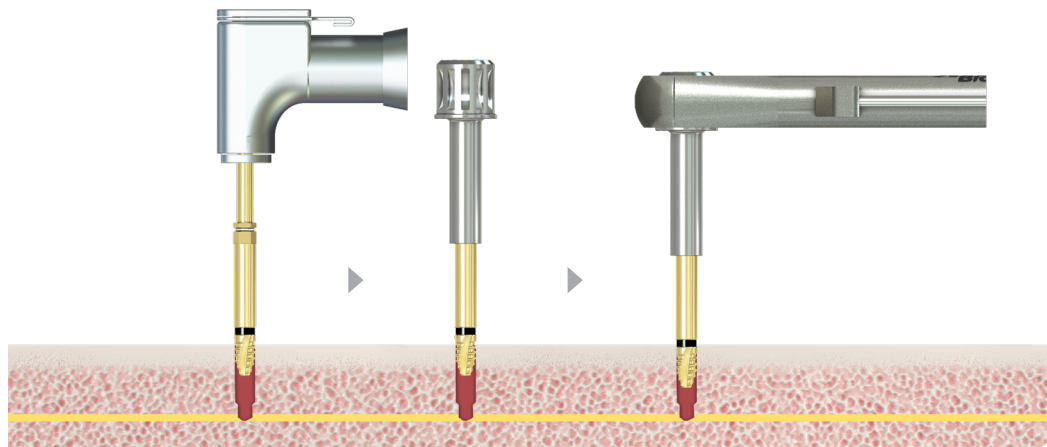

If the Cortical Tap Drill for contra angle gets stuck in a bone, continue tapping by connecting Ratchet Connector and Torque Ratchet



Neobiotech USA
2700 E Foothill Blvd. #121 Pasadena CA 91107
www.neobiotechusa.com 213-387-7704